



Shock Level Tester

SYSTEM COMPONENTS

USB Connection
9V Power Supply
Alligator Clips
Software and Manual

PART NUMBERS

➤ 2325-0358
SHOCK LEVEL TESTER

COMPUTER REQUIREMENTS

SDI offers high performance
Configure Computers that are
pre-installed with the Windows
operating system, USB Drivers
and applicable SDI software.
Windows 10 compatible
computer systems with one USB
port.

For more information on any of
our products or services please
visit us on the Web at:
www.sandiegoinstruments.com
Or contact us via email at
sales@sandiegoinstruments.com



FEATURES & BENEFITS

Displays peak value for both voltage
and amperage

Chart read out of Oscilloscope

Set duration test time up to 1 minute


Easily connects via USB

Alligator clips to secure an easy
connection

PRODUCT OVERVIEW

The SDI Shock Level Tester is a convenient solution to measure and verify that the levels of adverse stimulation are appropriate for the subject being tested. By programing your system to a desired current output, the Shock Level Tester Software should read the same to confirm accurate current output. The system will also display the waveform on the digital oscilloscope for analysis.

Outside Box Dimensions	7" x 6" x 2"
Control Unit Weight	2 lbs
Material Composition	ABS Plastic
Cable Clip Length	4 ft
Power Source	9V
Computer Connector	USB

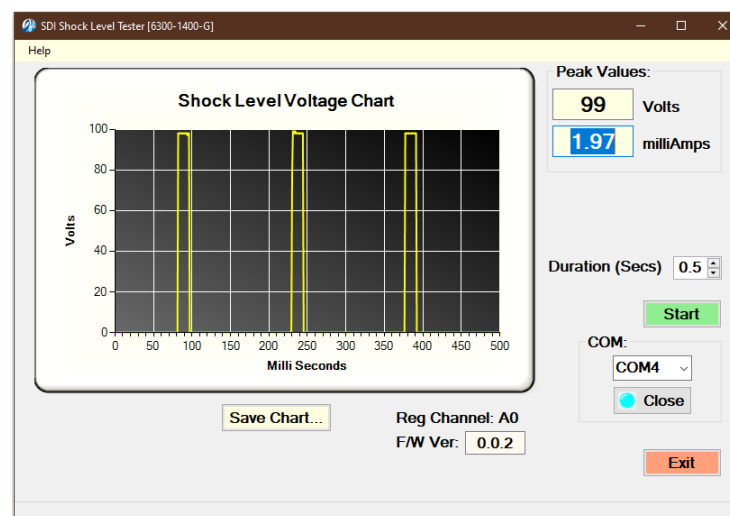

9155 Brown Deer Rd Suite 8
San Diego CA 92121
Phone 858-530-2600